

Allen.
Bow® Frame



Versatile – Use on your OR table, Jackson™ Table or Allen Flex Frame

The Allen Bow Frame is a platform used to perform less complex spine surgeries. The Bow Frame is more versatile than the existing Wilson™ Frame, as the single device attaches to standard OR tables, the Jackson Table and the Allen Flex Frame.

The Bow Frame also has a higher patient weight capacity than the Wilson Frame, of 500 lbs (227 kg). In addition, the Bow Frame Skin Care™ Cover disposables are comprised of moisture-absorbing foam, with a low shear backing.

Ideal for use with:

- Laminectomy
- Fusion procedures
- Microdisc procedures
- Decompressions
- Disc procedures
- Foraminotomy

At a glance

- Patient weight capacity: 500 lbs (227 kg)
- Fits the Jackson Table
- Fits the Allen Flex Frame
- Fits OR Tables
- Radiolucent
- Patient supports are adjustable to fit patients of various sizes.
- Patient supports adjust laterally up to 10 inches (25.4 cm) to relieve pressure on the patient's abdomen
- Lightweight for easy set-up
- Compact for easy storage
- Storage cart included for easy transport
- Skin Care Cover disposables are comprised of cosmetic foam, and a low shear backing designed to protect skin from friction and shear

About the Product

The Allen Bow Frame is used as a platform to perform spine surgeries. Its bow shape creates the desired reverse lordosis needed to open the intervertebral spaces. The frame is designed to fit standard OR Tables, and attach to the Allen Flex Frame, and Jackson™ Table. The Bow Frame has a patient weight limit of 500 lbs (227 kg).

The Bow Frame enables the patient's abdomen to hang free during procedures. It is also radiolucent and allows for unrestricted C-arm access. The removable crank is used to adjust the amount of flex and tension in the patient supports. Every Bow Frame comes with a storage cart.




Bow Frame



Crank Arm



Storage Cart

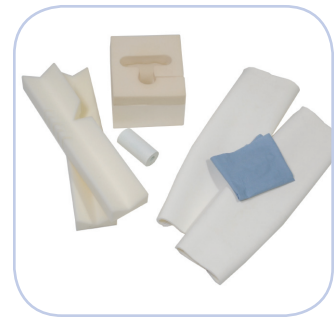

Bow Frame Skin Care Covers Kit

Technical Specifications

Patient Weight Capacity 500 lbs (227 kg)
Device Weight.....25 lbs (11.3 kg)

Allen Bow Frame Contents:

- Bow Frame
- Crank Arm
- Storage Cart
- Bow Frame Skin Care Covers (Includes 2 pad covers, 2 arm cradles, 1 foam block head positioner, 1 transverse bar cover, 1 base cover per set)

#A-70800 Allen Bow Frame

#A-70810 Skin Care Cover Disposables (Case of 6)

Allen® is a registered trademark of Allen Medical Systems.
Jackson™ and Wilson™ are trademarks of Mizuho Orthopaedic Systems.

Hill-Rom reserves the right to make changes without notice in design, specifications and models. The only warranty Hill-Rom makes is the express written warranty extended on the sale or rental of its products.

© 2012 Allen Medical, A Hill-Rom Company.
ALL RIGHTS RESERVED.

D-770575-A4



USA 800-433-5774
www.allenmedical.com